



PureRegen[®] Gel

Injectable Crosslinked Hyaluronic Acid Nasal Packing and Stent

PureRegen[®] Injectable Nasal Packing is formulated to be used as a nasal dressing and stent, separating nasal structures following nasal and sinus surgical procedures.



Controls Minimal Bleeding and Oozing

Clinically Tested Effective Adhesion Prevention

Degradable - No Need for Painful Traumatic Removal

Facilitates Wound Healing and Functional
Mucociliary Regeneration



800 998 8580
www.invotec.net

PureRegen® Gel

Injectable Crosslinked Hyaluronic Acid Nasal Packing and Stent

Comprehensive nasal wound healing and more efficient mucosal regeneration.

PureRegen® bioresorbable nasal packing and stent is a non-pyrogenic, transparent, specifically formulated viscoelastic gel of self-cross-linked hyaluronan derived from nonanimal sources. This patented formulation ensures the material maintains the best possible stability during its storage and clinical use. Unlike esterified HA, PureRegen® Gel has not been reported as osteogenic in the nasal cavities.

PureRegen® Gel absorbs body fluids and blood resulting in its expansion gradually up to 50% by volume. The volume of the injected gel secures its initial stability within sinus cavity and controls minimal bleeding by tamponade effect.


Predictable retention, maximum tissue coverage and adhesion prevention

Supports most favorable tissue separation


Patient comfort

Eliminates need for traumatic removal

Easy to use pre-filled syringe with malleable cannula

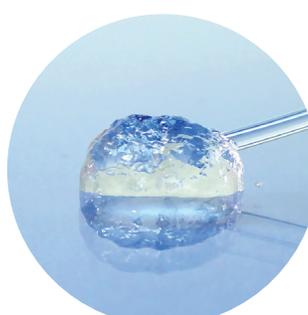

PureRegen® Gel custom formulated viscosity provides maximum tissue coverage, and predictable retention resulting in the best possible adhesion prevention.

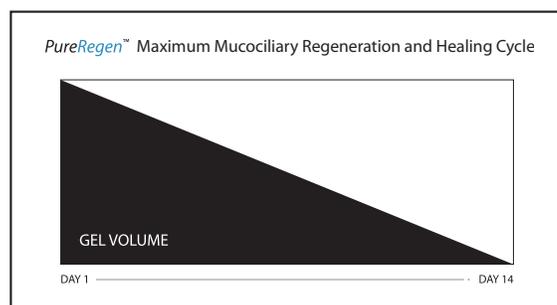

PureRegen® Gel jackets the nasal mucosa during the full 14-day post op inflammatory response period facilitating comprehensive healing and functional mucociliary regeneration.

Ordering Information

Catalog No.	Description	Quantity
PR-00005	PureRegen® Gel - 5 ml prefilled glass syringe with cannula	1 package



6833 PHILLIPS INDUSTRIAL BLVD.
JACKSONVILLE, FL 32256